

Application News



Statistical analysis of your database

ACTIVITY REPORT

Time frame
 Last months Date range from 12.10.2009 to 12.10.2010

Activities to be reported

TOTAL NUMBER OF STUDIES PEDIATRIC ECHO STUDIES ABDOMINAL AND SHARED SERVICE

ADULT ECHO STUDIES 0 - 2 years Abdominal

Cardiac studies 2 - 5 years Small parts

Contrast studies 5 - 10 years Ob/Gyn

Stress echo studies 10 - 15 years Contrast

2-level stress 15 - 20 years

3-level stress 20 years and older

4-level stress VASCULAR

5-level stress Carotid studies

With contrast Lower extremity vein

FETAL ECHO STUDIES Lower extremity arterial

0 - 12 weeks Upper extremity vein

12 weeks and older Upper extremity arterial

Fetal heart

Show activities for

Sonographer (operator)
 Lab

Template Standard (la

Activity Report

Reporting Period: From 13/10/2009 to 12/10/2010

Sonographer (operator)

	Total No Studies	Adult Echo Studies	Fetal Echo Studies	Pediatric Total
Total	50	48	0	2
ADM	9	7	-	2
DOK	1	1	-	-
GRL	4	4	-	-
ABC	14	14	-	-
LAN	3	3	-	-
SRV	6	6	-	-
Unknown	5	5	-	-
USR	5	5	-	-

Department

	Total No Studies	Adult Echo Studies	Fetal Echo Studies	Pediatric Total
Total	50	48	0	2
Echo	41	41	-	-
Cardiology	1	-	-	1
Development	2	2	-	-
Kardiologie	5	5	-	-
ULTRASOUND	1	-	-	1



Content

STATISTICAL ANALYSIS OF YOUR DATABASE	2
HOW TO PERFORM STATISTIC ABOUT YOUR EXAMS	4
Note	4
How to get there	4
Set up your selections	4
The time frame	4
Study type	5
How to get the data out	6
Preview and Send	6
Sent to printer	6
Sent to file	7

NOTE

This hand out is additional training material.
For more information please refer to the user manual and/or reference manual.



How to perform statistic about your exams

Note


This Tool is only available on EchoPAC version Bt09 or higher

How to get there


Go to **Config** (on the upper row)

Select the **Admin** Tab

Press on the **Audit** Tab

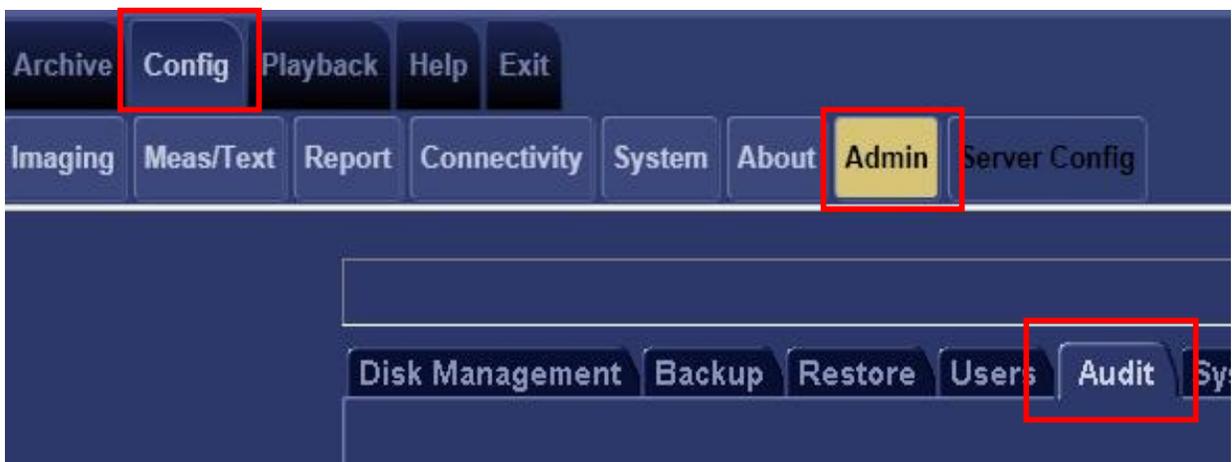


Select **View Activity** report

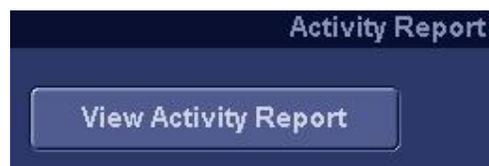

Set up your selections

The time frame

Select time frame you want to run the analysis through:
Either by using the selection of: days, weeks, months or years.

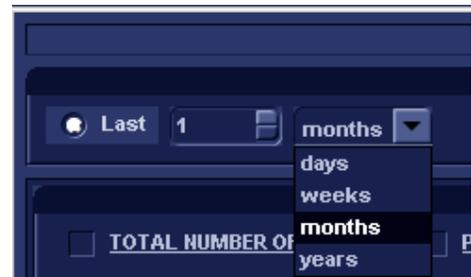

Or define a time frame with Date range from – to a specific date.


Study type

Select which kind of studies you'd like to include.
More than one selection can be done.

Separate the studies by selection of the following options

How to get the data out

Select in which format you would like to have your data.
Either as print out or stored as file.

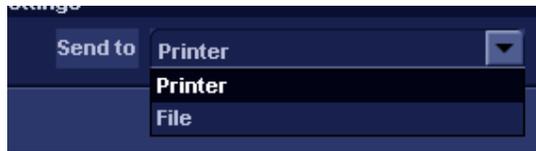

Preview and Send

If you would like to see the results, before sending them out, press the **Preview** button.



Activity Report

Reporting Period: From 12/10/2009 to 12/10/2010

Sonographer (operator)

	Adult Echo Studies	Total Cardiac Studies	Fetal Echo Studies	Pediatric Total	Vascular Total	Abdominal And Shared
Total	48	48	0	2	0	0
ADM	7	7	-	2	-	-
DOK	1	1	-	-	-	-
GRL	4	4	-	-	-	-
HFF	14	14	-	-	-	-
LAN	3	3	-	-	-	-
SRV	6	6	-	-	-	-
Unknown	5	5	-	-	-	-
USR	5	5	-	-	-	-

In order to start the print/store file press the **Send** button.



Sent to printer

When selected to send the analysis to the printer, press Send and the system will print the report on the printer that is setup for use in your EchoPAC.


Sent to file

If you selected to store the results in a file, press Send and a new window open up.

Select where to store the file.

You may change the file name, if needed.

Select the type of file you wanted to store: CHM (open with internet browser), PDF, and TXT

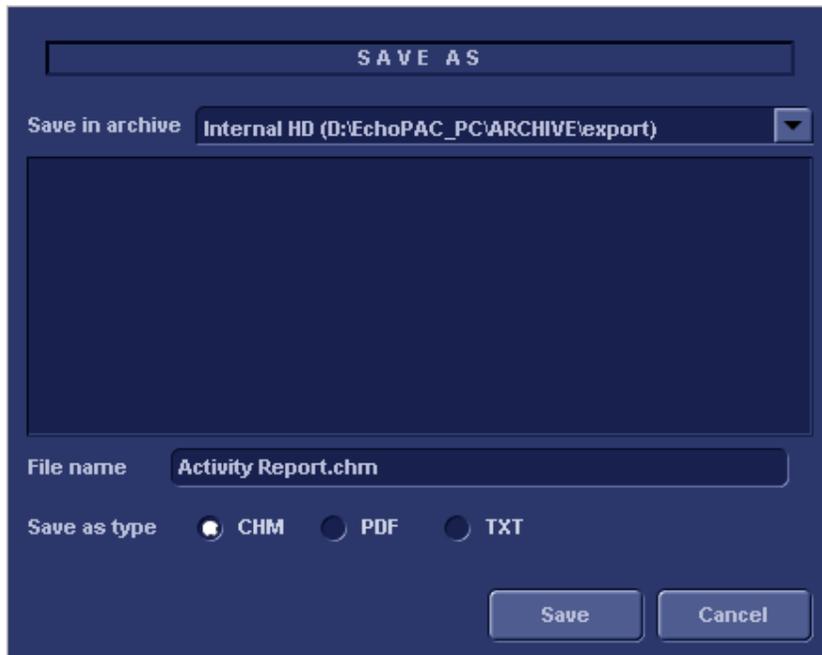

Close

Once you are finished press **Close** in order to end the session.